

## IRRIGATION LINES

### **Hose set for mechanically driven irrigation system**

Omnia offers a wide range of mechanical irrigation lines, specifically developed for the majority of Physiodispensers available on the market. Irrigation systems are supplied with a perforator, a roller for flow control, fittings, and an integrated peristaltic pump section.

Omnia mechanical irrigation systems offer the following options:


- a. single internal irrigation
- b. single external irrigation
- c. double internal and external irrigation

Each set for mechanical irrigation is packed in a sterile medical blister.

Each pack contains 10 pieces and each set is single use.

### **FEATURES**

- Perforator with air filter
- Roller to adjust the flow
- Pump section in medical grade silicone
- Total length of the line, from peristaltic pump to handpiece, 230 cm
- End tip in soft PVC

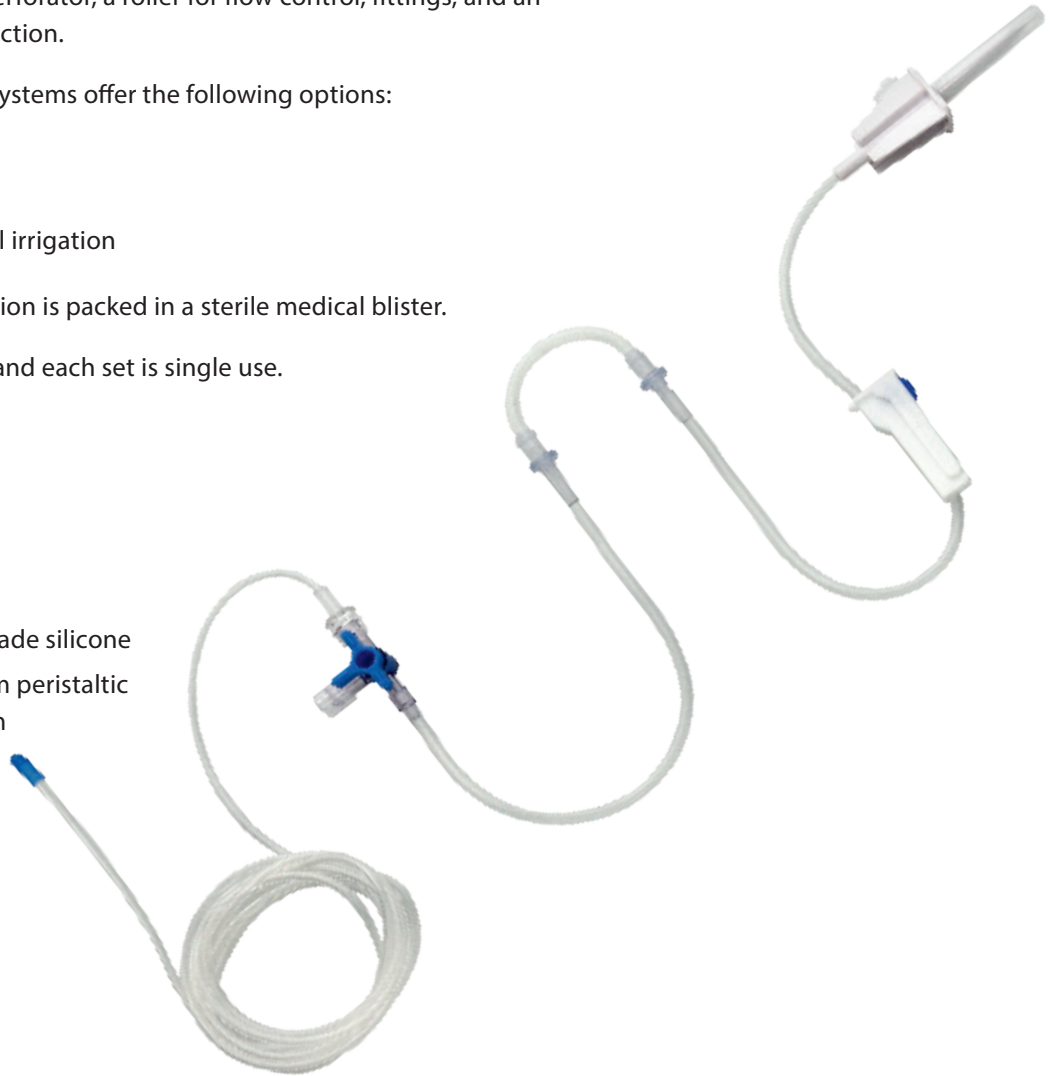




Disposable Medical Devices

## IRRIGATION LINES



<b>OMNIA® Irrigation Line Set For NSK Based Drilling Units</b> <b>SKU# OM32F0134</b>	<b>OMNIA® Irrigation Line Set For W&amp;H Based Drilling Units</b> <b>SKU# OM32F0139</b>
Irrigation set	Irrigation set
1/2Y	 1/2Y
Suitable for: NSK Surgic XT, Surgic XT Plus, Osseocision®, VarioSurg, Kyocera Implantor Neo, Aseptico Aeu 6000 series, Aseptico Aeu 7000 series, Bonart Piezo surgical system	Suitable for: Nobel Biocare80/ W&H Implantmed/ Implant Med new design 2005/ Implantmed SI 923/ Astratech Implant SI 923/ Surgical Motor System/Ace/Strauman/InnovaCorp. 3I/Centerpulse/Lifecore/ DU 900/FRIOS® UNIT S/ FRIOS® UNIT S/I, Drilltech/Osseoset 100/Osseoset 200/Osseocare, Dentium iCTmotor
10 pieces/ box	10 pieces/ box

1/2Y = internal and external irrigation with Y joint



= Tap for the connection of a second irrigation line to a second contrangle or for the connection of Omnirrigator